


VITAPM®9

Product information

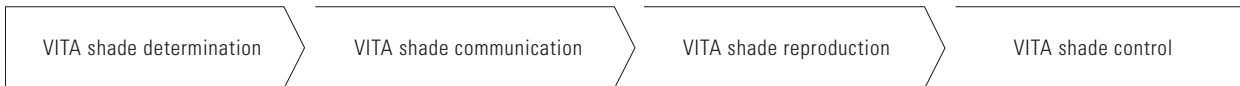

Date of issue: 06.15

VITA shade, VITA made.

VITA

VITAPM®9 Press Ceramic System

VITAVM®9 ADD-ON Individualization System

VITAPM[®]9 – guarantees economically efficient processing



VITAPM[®]9: All-in-one press ceramic for greater economic efficiency

VITA PM 9 (Pressable Material) is an all-in-one press ceramic and based on the fine-structure veneering ceramic VITA VM 9 which has proved its reliability in millions of restorations. VITA PM 9 covers a wide range of indications and ensures economic efficiency. It can be used for the substructure-free pressing technique (inlay/onlay, veneer and anterior crowns) as well as for pressing over zirconia substructures. The compact range of shades includes three degrees of translucency (O, T, HT) with 10 shade variations each.

The advantages

- Economically efficient fabrication: the wide range of indications (substructure-free pressing technique and overpressing technique) of VITA PM 9 enables economical and reliable fabrication of different types of restorations with the same pressing program, while at the same time maintaining efficient inventory levels.
- Clinical reliability: the VITA PM 9 ceramic was developed on the basis of the proven fine-structure veneering ceramic VITA VM 9 and ensures unsurpassed clinical reliability.
- Natural play of colors: in particular the highly translucent VITA PM 9 HT pellets demonstrate a convincing chameleon effect and guarantee a natural shade effect.
- Characteristics typical of natural teeth: VITA PM 9 ceramic demonstrates the characteristics of natural teeth, such as homogeneous surfaces, high degree of plaque resistance and enamel-like material properties.

The range of shades

The compact range of shades with 10 shade variations, each available in three different degrees of translucency (O = opaque, T = translucent, HT = high translucent), enables the fabrication of highly esthetic restorations.

VITA PM 9 is matched to the shades of the VITA SYSTEM 3D-MASTER shade system.

Overview of shades of O/T pellets:



	Shade variations
O pellets	0M1P, 0M2P, 1M1P, 1M2P, 2M1P, 2M2P 2M3P, 3M1P, 3M2P und 3M3P
T pellets	0M1P, 0M2P, 1M1P, 1M2P, 2M1P, 2M2P, 2M3P, 3M1P, 3M2P und 3M3P

Shade overview - T pellets::

	Shade variations
Dentin shades	0M2P, 1M1P, 1M2P, 2M2P, 3M2P
Enamel shades	EN0P, EN1P, EN2P, ENLP, ENDP

System compatibility

VITA PM 9 can be processed with all commercially available investment systems and ceramic press furnaces.

VITAVM[®]9 ADD-ON – guarantees esthetics of natural teeth

VITAVM[®]9 ADD-ON : Individualization materials for premium esthetics

Low-melting VITA VM 9 ADD-ON materials are now available especially for individualizing (layering technique) substructure-free VITA PM 9 restorations. These materials were developed on the basis of the proven VITA VM 9 veneering ceramic. They are available in eight different shade nuances and enable the dental technician to reproduce individual shade nuances quickly and reliably and obtain premium esthetics.



The advantages

- Natural shade effect: thanks to their outstanding translucency, enamel-like light refraction properties and harmonious combination with the VITA PM 9 base structure, VITA VM 9 ADD-ON materials ensure a natural play of colors.
- Reliable bonding: VITA VM 9 ADD-ON materials are optimally matched to the VITA PM 9 ceramic and form an excellent, distortion free sintered bond with the VITA VM 9 base structure. This ensures long-term load capacity.
- Economic efficiency: excellent stability during layering, minimized shrinkage behavior and first-class firing stability ensure efficient processing.
- Simple processing: VITA VM 9 ADD-ON materials produce a homogeneous material structure and thus guarantee good milling and polishing properties.



The range of shades

Enamel/ transparent materials	ADD1 (transparent), ADD2 (enamel light), ADD3 (enamel dark), ADD4 (whitish transparent)
Translucent materials	ADD5 (yellowish translucent), ADD6 (orange translucent), ADD7 (red translucent), ADD8 (blue translucent)

System compatibility

VITA VM 9 ADD-ON materials were especially developed for individualizing substructure-free VITA PM 9 restorations. The low-melting materials can be fired at 800 °C according to the working instructions with all commercially available ceramic furnaces.


VITA AKZENT Plus stains: – Optimized. Simple. Comprehensive.

With 19 shades, VITA AKZENT Plus offers dental practices and laboratories a complete range of products for numerous shade characterization options. VITA AKZENT Plus can not only be used for characterizing VITA PM 9 restorations; the stains are also suitable for all other dental ceramic materials, regardless of the restoration's CTE. These new fluorescent stains allow internal staining of restorations during layering as well as staining and glazing of external surfaces.

Depending on the user's preferred method of processing and the relevant area of application, VITA AKZENT Plus stains are available as powders and ready-to-use pastes. The glazing Body Stains and Glaze materials are also available as sprays.



With the unique VITA SYSTEM 3D-MASTER, all natural tooth shades can be systematically determined and perfectly reproduced.

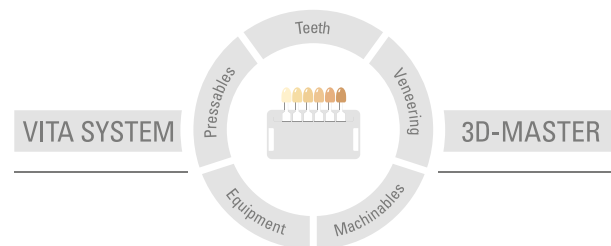

Please note: Our products must be used in accordance with the instructions for use. We accept no liability for any damage resulting from incorrect handling or usage. The user is furthermore obliged to check the product before use with regard to its suitability for the intended area of application. We cannot accept any liability if the product is used in conjunction with materials and equipment from other manufacturers that are not compatible or not authorized for use with our product. Furthermore, our liability for the accuracy of this information is independent of the legal basis and, in as far as legally permissible, shall always be limited to the value as invoiced of the goods supplied, excluding value-added tax. In particular, as far as legally permissible, we do not assume any liability for loss of earnings, indirect damages, ensuing damages or for third-party claims against the purchaser. Claims for damages based on fault liability (culpa in contrahendo, breach of contract, unlawful acts, etc.) can only be made in the case of intent or gross negligence. The VITA Modulbox is not necessarily a component of the product.

Date of issue of this information: 06.15

After the publication of these information for use any previous versions become obsolete. The current version can be found at www.vita-zahnfabrik.com

VITA Zahnfabrik has been certified in accordance to the Medical Device Directive and the following products bear the CE mark  0124 :

VITAVM[®]9 · VITAPM[®]9 · VITA AKZENT[®] Plus

VITA

VITA Zahnfabrik H. Rauter GmbH & Co.KG
Spitalgasse 3 · D-79713 Bad Säckingen · Germany
Tel. +49(0)7761/562-0 · Fax +49(0)7761/562-299
Hotline: Tel. +49(0)7761/562-222 · Fax +49(0)7761/562-446
www.vita-zahnfabrik.com · info@vita-zahnfabrik.com
 facebook.com/vita.zahnfabrik